

Multi-functional digital meridian physiotherapy instrument

■ Dredge the meridians ■ Relieve fatigue ■ Harmonize Qi and Blood ■ Physical fitness

description



The key to your energy.

Multi-functional digital meridian physiotherapy instrument DIGITAL PHYSIOTHERAPY MACHINE

Medium frequency physiotherapy is a kind of treatment method "from external treatment". It works through meridians, acupoints. Use certain techniques for physiotherapy to pass through the meridians, regulate the blood, and make the yin and Yang relatively balanced. The function of the viscera tends to be reconciled, so as to achieve the purpose of strengthening physical fitness.

Product introduction

Research and development center based on the principle of traditional Chinese medicine, using international advanced microcomputer electronic technology, research The utility model relates to a multifunctional microbrain electronic physiotherapy product with high technology content. It passes through the micro Computer controlled low wave simulation of Chinese medicine acupuncture, scraping, cupping, massage techniques, Through the human body meridian system conditioning the whole body, play through the meridian and Qi and blood, coordination Yin Jin, strong body The health function of the body.


Dredge meridian and Qi and blood to relieve fatigue and physical fitness

Shujing Huoluo Tongli joint Huoxue blood stasis wind scattered knot

Xingqi Huazheng Fuzheng dispel evil spirits to relieve fatigue and strengthen physical fitness

Product characteristics

- *Hammer, massage, massage, kneading, slapping, scraping, fingering, automatic mode, A variety of massage instruments in one
- *Can be applied to shoulders, waist, knees, hands and feet. Multiple parts of the leg, etc..
- *Easy to use, easy to operate, easy to use and easy to deliver
- *Pure natural therapy with no side effects or safety hazards

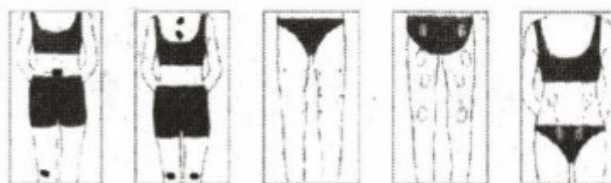

Multi-functional digital meridian physiotherapy instrument DIGITAL PHYSIOTHERAPY MACHINE

Product construction



host



Output hole Power DC hole Output hole



Physiotherapy shoes



A physiotherapy patch



AAA Batteries



Electric source



USB Wires



wire

Warm tips.

For the sake of your health and that of your family, please consult your doctor.

Patients with malignant tumors

Pregnant or menstruating women

Heart damage.

People with acute illnesses,

Patients undergoing medical treatment, especially those who feel abnormal body

Multi-functional digital meridian physiotherapy instrument DIGITAL PHYSIOTHERAPY MACHINE

method of application

1. Open the battery cover, follow the direction indicated, correctly insert the No. 3 5 battery and close the battery cover
2. Please insert the wire plug into the output socket L of the physiotherapy instrument, and the other end is correct with the physiotherapy slippers or physiotherapy patches. Connection.
3. Please put your feet on physical therapy shoes at the same time or paste the patch on the area where physiotherapy is needed(the part that requires physiotherapy) Please clean it first).
4. Long Press Open Key(upper left side of product) Open Massager
5. This machine is equipped with hammer, massage, massage, kneading, tapping, scraping, fingering, and eight physiotherapy modes in automatic mode. Free to choose according to individual needs.
6. This machine is equipped with physical therapy for shoulders, knees, hands and feet, legs, etc.. It can be used according to personal needs. By choice.
7. Strength + is the massage strength control key, which can be regulated by the country according to individual needs, if this machine No load detected, 5 SEC turn off. When the speed button is attached, the Si stupid D bit on the LCD display Sets the indicator for the intensity file position of 1-15.
8. When the product is working, please paste both sides of the physiotherapy label firmly at the same time to prevent detection of light or no load Load Hats +
9. Press the key machine to stop working.
10. This machine has a 5-30 minute automatic shutdown function. The LED screen "time" randomly displays the remaining working hours when the LED screen is used. The screen shows "0 time" and the machine automatically stops working.

Fault Analysis and Processing

Exceptional situation	reason	processing
Power's down.	Is the socket plug in good contact?	Check power connection plug and power switch. Is it normal?
		Is the battery in good contact with the cells
Automatic shutdown	Physical therapy patches fall off or Poor connection contact	Please clean the physiotherapy area and put the patch on.Reload Check the power cable.
Physical therapy is weak.	Operation with or without voltage instability	Turn off the power. Turn it back on.
Some parts feel uncomfortable	Whether use time is too long	Over 30 minutes per use
	Prolonged physiotherapy Same area.	Please change the other parts.

Multi-functional digital meridian physiotherapy instrument DIGITAL PHYSIOTHERAPY MACHINE

Attention

- * Before use, please paste the physiotherapy patch on the skin before turning on the power supply.
- * During use, do not move the patch up or onto other people.
If you feel any discomfort during use, stop using it immediately.
- * Do not use power sources other than rated voltage
- * Please do not disassemble and repair this machine except for the repair technicians appointed by the company Managed and modified.
- * Do not use in the bathroom, kitchen and other humid or hot places.

Custody method

- * If there is dirt or dust, please use a soft cloth with neutral detergent to wipe and air Dry naturally.
- * If the machine stops working too long, please remove the battery and save it separately.
- * When not in use, do not store in hot places such as stove or kitchen
- * Do not store in wet places when not in use



warning

Do not use tiana water, sulfuric acid and other corrosive liquid wipe.

Multi-functional digital meridian physiotherapy instrument DIGITAL PHYSIOTHERAPY MACHINE

After-Sales service and producer responsibility

- *The warranty period of this product is one year, calculated from the date of purchase.
- *During the quality assurance period, if non-human damage can be sent to the purchase point with invoice and warranty card to enjoy Free Maintenance Services
- *Subject to the following conditions, materials will be collected as appropriate during the valid warranty period Principal or service charge
 - Failure to show purchase invoice and warranty card
 - Damage caused by improper or incorrect use
 - Failure due to self-repair or modification
 - Damage caused by irresistible factors

Technical parameters

Product Name	Multi-functional medium frequency physiotherapist		
Rated voltage	4.5V	Rated working hours	5-30minute
Rated power	0.15W	weight	176G
Accretion	18cm*12.5cm*5cm	standard of execution	GB4706.1-1998
			GB4706.10-2003



This product has been rigorously tested and proved to be safe and reliable

To ensure that you can use this product correctly, please read this specification carefully and be careful
Keep this manual for future queries.

We reserve the right to make improvements to our products and specifications without prior notice.

The appearance of the product is based on the physical object.

